



## SynergEyes® A, SynergEyes M, SynergEyes KC, SynergEyes ClearKone® SynergEyes PS

### HDS (paflufocon D hem-iberfilcon A) HYBRID CONTACT LENSES FOR DAILY WEAR



**CAUTION: FEDERAL LAW RESTRICTS THIS DEVICE TO SALE BY OR ON THE ORDER OF A LICENSED EYE CARE PRACTITIONER.**

*Important: Please read carefully and keep this information for future use. This package insert is intended for the eye care practitioner but should be made available to patients upon request. The eye care practitioner should provide the patient with patient instructions that pertain to the patient's prescribed lens.*

#### SynergEyes A and SynergEyes M

##### **DESCRIPTION:**

SynergEyes A and SynergEyes M HDS (paflufocon D hem-iberfilcon A) Hybrid Contact Lenses for hyperopia and myopia with and without astigmatism, and presbyopia, are for daily wear.

<b>SynergEyes A and SynergEyes M Lens Parameters Available</b>	<b>RGP Center</b>	<b>Soft Skirt</b>	<b>Overall Lens</b>
Diameter (Approximate)	-	-	14.5 mm
Base Curve Radius	7.0 to 8.5 mm	-	-
Skirt Curve	-	0.7 to 1.3 mm flatter than BC	-
Posterior Optical Zone Diameter	7.80 mm	-	-
Center Thickness Range (for EP add 0.09 mm)	0.16 to 0.34 mm	-	-
Dioptric Power	+20.00 to -20.00 D	-	-
Add Power – Multifocal Only	+0.75 to +4.00 D	-	-
Zone – Multifocal Only	1.0 to 5.0 mm	-	-
<b>Lens Physical/Optical Properties</b>			
Refractive Index	1.442 Nd at 25°C	-	-
Spectral Transmittance (D&C Green 6) (380 to 780 nm) *	>90%	-	-
Wetting Angle (Receding Angle)	42°	-	-
Specific Gravity	1.10	-	-
Oxygen Permeability†	100	9.3	-
Water Content	<1%	27%	-

\* Ophthalmic Optics - Contact Lenses - Part 3: Measurement Methods, ISO 18369-3

† Ophthalmic Optics - Contact Lenses - Part 4: Physicochemical Properties of Contact Lens Materials, ISO 18369-4

##### **INDICATIONS FOR USE:**

SynergEyes A (paflufocon D hem-iberfilcon A) and SynergEyes M (paflufocon D hem-iberfilcon A) Hybrid Contact Lenses are indicated for use in the correction of hyperopic, myopic, and astigmatic refractive error including presbyopia, in aphakic and not aphakic, non-diseased eyes. The lenses are indicated for the correction of up to +20.00 and -20.00 D in eyes with astigmatism up to 6.00 D. For presbyopia, add powers between +1.00 and +4.00 D. The lenses may be prescribed for daily wear with removal for cleaning and disinfection (chemical, not heat) prior to reinsertion as recommended by the eye care professional.

#### SynergEyes KC and SynergEyes ClearKone

SynergEyes KC (paflufocon D hem-iberfilcon A) and SynergEyes *ClearKone* (paflufocon D hem-iberfilcon A) Hybrid Contact Lenses provide refractive error correction for irregular astigmatism/keratoconus when worn for daily wear.

ClearKone only: Generally, the central base curve is chosen to be deeper than the curvature of the flattest central cornea meridian. The lens is fitted to allow this zone to clear the central corneal apex. The lens is intended to maintain alignment over the keratoconic anterior geometry. In the case of corneal astigmatism increased clearance will be observed in the peripheral portion of the steepest meridian and some bearing may appear in the paracentral region of the flattest meridian.

KC and ClearKone: The central rigid portion of the lens measures 8.4 mm. The transition to the peripheral spherical zone begins outside the rigid – soft junction in a seamless fashion. A posterior peripheral bevel is present and terminates at the lens edge. The lens diameter is held constant at 14.5 mm. Post lens tear exchange is facilitated by a pumping action upon blinking and assisted by the minimal edge clearance provided by the peripheral bevel. The maximum central thickness of minus power lenses and the maximum junction thickness of plus power lenses is 0.20 mm. The center thickness of minus power lenses and the anterior optic zone junction thickness of plus power lenses reduces as the lens power increases to hold a constant rigid – soft junction thickness, soft skirt thickness and lens edge thickness. The SynergEyes KC and ClearKone Hybrid contact lens edge is pre-specified and equivalent in all lenses regardless of their other parameters. The anterior central curve is selected to provide any necessary optical power to correct spherical refractive error not corrected by the optical effect of the posterior base curve and the tear lens formed between it and the cornea. As with rigid gas permeable lenses there may be residual astigmatism uncorrected by the lenses. The amount of residual astigmatism may be estimated by comparison of the corneal and refractive astigmatism. Eyes with near equal corneal and refractive astigmatism are not expected to demonstrate residual astigmatism. Eyes with a disparity between corneal and refractive astigmatism of greater than 0.75 D may demonstrate residual astigmatism.

<b>SynergEyes KC and SynergEyes ClearKone Lens Parameters Available</b>	<b>RGP Center</b>		<b>Soft Skirt</b>	<b>Overall Lens</b>
	<b>KC</b>	<b>ClearKone</b>	<b>KC &amp; ClearKone</b>	
Diameter (Approximate)	-	-	-	14.5 mm
Base Curve Radius	5.3 to 8.5 mm	-	-	-
Vault	-	0 to 750 µ	-	-
Skirt Curve	-	-	7.6 to 9.1 mm	-
Posterior Optical Zone Diameter	9.0 mm	6.0 mm	-	-
Center Thickness Range (for EP add 0.09 mm)	0.18 to 0.34 mm	0.21 to 0.39 mm	-	-
Dioptric Power	+20.00 to -20.00 D	+20.00 to -20.00 D	-	-
<b>Lens Physical/Optical Properties</b>				
Refractive Index	1.442 Nd at 25°C		-	-
Spectral Transmittance (D&C Green 6) (380 to 780 nm) *	>90%		-	-
Wetting Angle (Receding Angle)	42°		-	-
Specific Gravity	1.10		-	-
Oxygen Permeability†	100		9.3	-
Water Content	<1%		27.1% ± 2.5%	-

##### **INDICATIONS FOR USE:**

SynergEyes KC (paflufocon D hem-iberfilcon A) and SynergEyes ClearKone (paflufocon D hem-iberfilcon A) Hybrid Contact Lenses for keratoconus are indicated for use in the correction of eyes with refractive errors that include hyperopia and myopia that manifest irregular astigmatism, in aphakic and not aphakic, and otherwise non-diseased eyes. The lenses are indicated for the correction of up to +20.00 and -20.00 D in eyes with irregular astigmatism up to 6.00 D. The lenses may be prescribed for daily wear with removal for cleaning and disinfection (chemical, not heat) prior to reinsertion as recommended by the eye care professional.

#### SynergEyes PS

SynergEyes PS (paflufocon D-hem-iberfilcon A) Hybrid Contact Lenses for Post-Surgical Refractive Error and Trauma provide refractive error correction for post-surgical or traumatic corneal conditions and are intended for daily wear.

The central rigid portion of the lens measures 8.4 mm. The transition to the peripheral spherical zone begins outside the rigid – soft junction in a seamless fashion. A posterior peripheral bevel is present and terminates at the lens edge. The lens diameter is held constant at 14.5 mm. Post lens tear exchange is facilitated by a peristaltic pumping action upon blinking and assisted by the minimal edge clearance provided by the peripheral bevel. The maximum central thickness of minus power lenses and the maximum junction thickness of plus power lenses is 0.30 mm. The center thickness of minus power lenses and the anterior optic zone junction thickness of plus power lenses reduces as the lens power increases to hold a constant rigid – soft junction thickness, soft skirt thickness and lens edge thickness. The SynergEyes PS Hybrid contact lens edge is pre-specified and equivalent in all lenses regardless of their other parameters. The anterior central curve is selected to provide any necessary optical power to correct spherical and cylindrical refractive error not corrected by the optical and mechanical effect of the posterior base curve and the tear lens formed between it and the cornea. As with rigid gas permeable lenses there may be residual astigmatism uncorrected by the lenses. The amount of residual astigmatism may be estimated by comparison of the corneal and refractive astigmatism. Eyes with near equal corneal and refractive astigmatism are not expected to demonstrate residual astigmatism. Eyes with a disparity between corneal and refractive astigmatism of greater than 0.75 D may demonstrate residual astigmatism and may require the non-rotating design with a toric anterior surface.

<b>SynergEyes PS Lens Parameters Available</b>	<b>RGP Center</b>	<b>Soft Skirt</b>	<b>Overall Lens</b>
Diameter (Approximate)	-	-	14.5 mm
Base Curve Radius	7.1 to 9.0 mm	-	-
Skirt Curve	-	8.3 to 8.6 mm	-
Posterior Optical Zone Diameter	6.5 mm	-	-
Center Thickness Range (for EP add 0.09 mm)	0.16 to 0.34 mm	-	-
Dioptric Power	+20.00 to -20.00 D	-	-
<b>Lens Physical/Optical Properties</b>			
Refractive Index	1.442 Nd at 25°C		-
Spectral Transmittance (D&C Green 6) (380 to 780 nm) *	>90%		-
Wetting Angle (Receding Angle)	42°		-
Specific Gravity	1.10		-
Oxygen Permeability†	100	9.3	-
Water Content	<1%	27%	-

##### **INDICATIONS FOR USE:**

SynergEyes PS (paflufocon D-hem-iberfilcon A) Hybrid Contact Lenses are indicated for use in the correction of eyes with refractive errors resulting from corneal surgery or trauma including hyperopia and myopia, astigmatism and irregular astigmatism in aphakic and not aphakic, non-diseased eyes. The lenses are indicated for the correction of up to +20.00 and -20.00 D in eyes with astigmatism up to 6.00 D. The lenses may be prescribed for daily wear with removal for cleaning and disinfection (chemical, not heat) prior to reinsertion as recommended by the eye care professional.

### All HDS Hybrid Contact Lens Models

##### **DESCRIPTION:**

All SynergEyes HDS Hybrid Contact Lenses are manufactured from rigid gas permeable material (paflufocon D) and SynergEyes poly-hema material (hem-iberfilcon A). The lenses are manufactured from Paragon HDS® 100 with a poly-hema hydrogel skirt. The lens center provides the optics of rigid gas permeable lenses while the hydrogel skirt provides the stability and performance of hydrogel lenses. Greater attention must be directed toward fitting the lens than with soft hydrogel lenses. The base curve of the lens is modulated to provide an optimum central lens-cornea relationship, the skirt radius is modulated to provide an optimum scleral relationship, and the power of the lens is modulated to provide the desired refractive correction.

The center material is a thermoset fluorosilicone acrylate copolymer derived primarily from siloxane acrylate, trifluoroethyl methacrylate, and methylmethacrylate. The peripheral skirt material is a poly-hema hydrophilic copolymer. The lenses are available as lathe-cut contact lenses with a blue tint in the rigid central material. The blue material contains D&C Green No. 6. The poly-hema skirt is clear.

RGP Material:	The central aspherical or spherical zone
	The intermediate aspherical or spherical zone
Poly-HEMA Material:	The peripheral anterior edge taper and posterior bevel
	An edge terminus smoothly joining the anterior taper to the posterior bevel

##### **ACTIONS:**

SynergEyes HDS Hybrid Contact Lenses act as a refracting medium to focus light rays on the retina.

##### **CONTRAINDICATIONS (REASONS NOT TO USE):**

###### **DO NOT USE HDS Hybrid Contact Lenses when any of the following conditions exist:**

- Acute and subacute inflammation or infection of the anterior chamber of the eye.
- Severe insufficiency of lacrimal secretion (dry eyes).
- If eyes become red or irritated.
- Corneal hypoesthesia (reduced corneal sensitivity).
- Any systemic disease that may affect the eye or be exaggerated by wearing contact lenses.
- Allergic reactions of ocular surfaces or adnexa that may be induced or exaggerated by wearing contact lenses or use of contact lens solutions.
- Allergy to any ingredient, such as mercury or Thimerosal, in a solution which is to be used to care for HDS Hybrid Contact Lenses.
- Any active corneal infection (bacterial, fungal, or viral).
- **SynergEyes A and SynergEyes M:** Any eye disease, injury, or abnormality that affects the cornea, conjunctiva, or eyelids.
- **SynergEyes KC and SynergEyes ClearKone:** Any eye disease (excluding Keratoconus), injury, or abnormality that affects the cornea, conjunctiva, or eyelids.
- **SynergEyes PS:** Any eye disease or abnormality that affects the cornea, conjunctiva, or eyelids except corneal trauma deemed to benefit from the device.

##### **WARNINGS:**

**EYE PROBLEMS, INCLUDING CORNEAL ULCERS, CAN DEVELOP RAPIDLY AND LEAD TO LOSS OF VISION. THE PATIENT SHOULD BE INSTRUCTED TO IMMEDIATELY REMOVE THE LENSES, AND PROMPTLY CONTACT THEIR EYE CARE PRACTITIONER IF THEY EXPERIENCE:**

- Eye Discomfort
- Excessive Tearing
- Vision Changes
- Loss of Vision
- Eye Redness
- Or Other Eye Problem

**Patients should be advised of the following warnings pertaining to contact lens wear:**

- **Problems with contact lenses and lens care products could result in serious injury to the eye.** It is essential that patients follow the eye care practitioner's directions and all labeling instructions for proper use of their lenses and lens care products, including the lens case. Eye problems, including corneal ulcers, can develop rapidly and lead to loss of vision.
- Daily wear lenses are not indicated for overnight wear, and patients should be instructed not to wear lenses while sleeping. Clinical studies have shown that the risk of serious adverse reactions is increased when these lenses are worn overnight.
- Studies have shown that contact lens wearers who are smokers have a higher incidence of adverse reactions than nonsmokers.
- If a patient experiences eye discomfort, excessive tearing, vision changes, loss of vision, or redness of the eye, the patient should be instructed to immediately remove the lens and promptly contact their eye care practitioner.
- All contact lens wearers must see their eye care practitioner as directed.

**Patients should be advised of the following instructions for use and warnings pertaining to contact lens wear:**

<b>Soaking and Storing Lenses:</b>	
Instruction for Use	<ul style="list-style-type: none"> <li>• Patients should be instructed to use only approved, fresh contact lens disinfecting solution each time they soak (store) their lenses.</li> </ul>
Warning	<ul style="list-style-type: none"> <li>• Patients should be instructed to not reuse or “top off” old solution left in their lens case since solution reuse reduces effective lens disinfection and could lead to severe infection, vision loss, or blindness. “Topping-Off” is the addition of fresh solution to solution that has been sitting the case.</li> </ul>
<b>Rub and Rinse Time:</b>	
Instruction for Use	<ul style="list-style-type: none"> <li>• Patients should be instructed to rub and rinse their lenses according to the recommended lens rubbing and rinsing times in the labeling of their disinfecting solution to adequately disinfect their lenses.</li> </ul>
Warning	<ul style="list-style-type: none"> <li>• Patients should be instructed to rub and rinse their lenses for the recommended amount of time to help prevent serious eye infections.</li> <li>• Patients should be instructed to never use water, saline solution, or rewetting drops to disinfect their lenses. These solutions will not disinfect their lenses. Not using the recommended disinfectant can lead to severe infection, vision loss, or blindness.</li> </ul>
<b>Lens Case Care:</b>	
Instruction for Use	<ul style="list-style-type: none"> <li>• Patients should be instructed to empty and clean contact lens cases with digital rubbing using fresh, sterile disinfecting solutions/contact lens cleaner. Never use water. It is recommended that cleaning is followed by rinsing with fresh, sterile disinfecting solutions (never use water) and wiping the lens cases with a fresh, clean tissue. After use, never air-dry or recap the lens case lids without performing proper cleaning methods. If air drying, be sure that no residual solution remains in the case before allowing it to air dry.</li> <li>• Patients should be instructed to replace their lens case according to the directions given by the eye care practitioner or the labeling that came with the case.</li> </ul>
Warning	<ul style="list-style-type: none"> <li>• Contact lens cases can be a source of bacterial growth.</li> <li>• Patients should be instructed to not store lenses or rinse the lens case with water or any non-sterile solution. Patients should be instructed to only use fresh disinfecting solution so they do not contaminate their lenses or lens case. Use of non-sterile solution can lead to severe infection, vision loss, or blindness.</li> </ul>
<b>Water Activity:</b>	
Instruction for Use	<ul style="list-style-type: none"> <li>• Patients should be instructed to not expose their contact lenses to water while wearing them.</li> </ul>
Warning	<ul style="list-style-type: none"> <li>• Water can harbor microorganisms that can lead to severe infection, vision loss, or blindness. Patients should be instructed that if their lenses have been submersed in water when swimming in pools, lakes, or oceans, they should discard them and replace them with a new pair. Patients should be instructed to ask the eye care practitioner for recommendations about wearing their lenses during any activity involving water.</li> </ul>
<b>Discard Date on Contact Lens Disinfecting Solution Bottle:</b>	
Instruction for Use	<ul style="list-style-type: none"> <li>• Patients should be instructed to discard any remaining solution after the recommended time period indicated on the bottle of disinfecting solution used for disinfecting and soaking their contact lenses.</li> <li>• Patients should be instructed that the Discard Date refers to the time they can safely use contact lens care product after the bottle has been opened; it is not the same as the expiration date, which is the last date that the product is still effective before it is opened.</li> </ul>
Warning	<ul style="list-style-type: none"> <li>• Patients should be instructed that using their disinfecting solution beyond the discard date could result in contamination of the solution and can lead to severe infection, vision loss, or blindness.</li> <li>• To avoid contamination, patients should be instructed: DO NOT touch the tip of the container to any surface. Replace cap after using.</li> <li>• To avoid contaminating their solution, patients should be instructed: DO NOT transfer to other bottles or containers.</li> </ul>

##### **PRECAUTIONS:**

###### **Special Precautions for Eye Care Practitioners:**

- Clinical studies demonstrated that contact lenses manufactured from the SynergEyes material are safe and effective for daily wear. Due to the small number of patients enrolled in clinical investigation of lenses, all refractive powers, design configurations, or lens parameters available in the HDS Hybrid Contact Lens material were not evaluated in significant numbers. Consequently, when selecting an appropriate lens design and parameters, the eye care practitioner should consider all characteristics of the lens that can affect lens performance and ocular health, including oxygen permeability, wettability, central and peripheral thickness, and optic zone diameter. The potential impact of these factors on the patient's ocular health should be carefully weighed against the patient's need for refractive correction; therefore, the continuing ocular health of the patient and lens performance on the eye should be carefully monitored by the prescribing eye care practitioner.
- Patients who wear aspheric contact lenses to correct presbyopia may not achieve the best corrected visual acuity for either far or near vision. Visual requirements vary with the individual and should be considered when selecting the most appropriate type of lens for each patient.
- Aphakic patients should not be fitted with HDS Hybrid Contact Lenses until the determination is made that the eye has healed completely.
- When using fluorescein in the fitting evaluation of Hybrid Contact Lenses; either high or standard molecular weight fluorescein can be used. The soft skirt of the lens may absorb this dye and become discolored. Whenever fluorescein is used in eyes, the eyes should be flushed with a sterile saline solution that is recommended for in-eye use. Lens discoloration may be reversed following the ISO guidelines for hygienic management of multipatient use trial lenses; if dye is still present, increase the hydrogen peroxide soak time.
- Before leaving the eye care practitioner's office, the patient should be able to promptly remove lenses or should have someone else available who can remove the lenses for him or her.
- Eye care practitioners should instruct the patient to remove the lenses immediately if the eyes become red or irritated.
- The patient should be instructed to always discard lenses after the recommended wearing schedule for a period of up to six (6) months or as recommended by the eye care practitioner.

**Eye care practitioners should carefully instruct patients about the following care regimen and safety precautions:**

- Different solutions cannot always be used together, and not all solutions are safe for use with all lenses. Use only recommended solutions. Chemical disinfection solutions should not be used with heat.
- Always use fresh, unexpired lens care solutions. Do not change lens care solutions without consulting the eye care practitioner.
- Always follow directions in the package inserts for the use of contact lens solutions.
- Do not heat the wetting/soaking solution and lenses. Keep away from extreme heat.
- Use only a chemical (NOT HEAT) lens care system. Use of a heat (thermal) care system can warp the center of the HDS Hybrid Contact Lenses.
- Sterile unpreserved solutions, when used, should be discarded after the time specified in the labeling directions.

- Do not use saliva, tap water, or anything other than recommended solutions for lubricating or wetting lenses.
- Always keep the lenses completely immersed in the recommended storage solution when the lenses are not being worn. Prolonged periods of drying may damage the poly-hema lens skirt. If the lens skirt becomes dried out, follow the Care for a *Dehydrated Lens* section of this instruction sheet.
- If the lens sticks (stops moving) on the eye, follow the *Care for a Sticking Lens* section of this instruction sheet. The lens should move freely on the eye for the continued health of the eye. If non-movement of the lens continues, the patient should be instructed to IMMEDIATELY consult his or her eye care practitioner.
- Always wash and rinse hands before handling lenses. Do not get cosmetics, lotions, soaps, creams, deodorants, or sprays in the eyes or on the lenses. It is best to put on lenses before putting on makeup. Water-based cosmetics are less likely to damage lenses than oil-based products.
- Do not touch contact lenses with fingers or hands if hands are not free of foreign materials; microscopic scratches of the lenses may occur, causing distorted vision and/or injury to the eye.
- Carefully follow handling, insertion, removal, cleaning, disinfecting, storing, and wearing instructions in the Patient Instruction Booklet for the HDS Hybrid Contact Lenses and those prescribed by the eye care practitioner. The Patient Instruction Booklet is available at SynergEyes.com.
- Never wear lenses beyond the period recommended by the eye care practitioner.
- If aerosol products, such as hairspray, are used while wearing lenses, exercise caution and keep eyes closed until the spray has settled.
- Avoid all harmful or irritating vapors and fumes while wearing lenses.
- Always handle lenses carefully and avoid dropping them. Never use tweezers or other tools to remove lenses from the lens container unless specifically indicated for that use. Pour the lens into the hand. Do not touch the lens with fingernails.
- Ask the eye care practitioner about wearing lenses during sporting activities.
- Inform the doctor (health care practitioner) about being a contact lens wearer.
- Always inform the employer about being a contact lens wearer; some jobs may require use of eye protection equipment or may require that the patient not wear contact lenses.
- Always discard lenses worn on a frequent replacement schedule after the recommended wearing schedule prescribed by the eye care practitioner.
- Always contact the eye care practitioner before using any medicine in the eyes.
- Contact lenses are medical devices and parameters or lens type should not be changed without consulting an eye care practitioner.
- Contact lenses should never be shared between users. Sharing contact lenses with someone else increases the risks for serious eye conditions.
- Never wear lenses with bubbles present.
- As with any contact lens, follow-up visits are necessary to assure the continuing health of the patient’s eyes. The patient should be instructed as to a recommended follow-up schedule.

**ADVERSE REACTIONS:**

**Patients should be informed that the following problems may occur:**

- Eyes stinging, burning, itching (irritation), or other eye pain
- Comfort is less than when lens was first placed on eye
- Feeling of something in the eye such as a foreign body or scratched area
- Excessive watering (tearing) of the eyes
- Unusual eye secretions
- Redness of the eyes
- Reduced sharpness of vision (poor visual acuity)
- Blurred vision, rainbows, or halos around objects
- Sensitivity to light (photophobia)
- Dry eyes

If the patient notices any of the above, he or she should be instructed to IMMEDIATELY REMOVE THEIR LENSES.

- If the discomfort or problem stops, then look closely at the lens; IF THE LENS IS IN ANY WAY DAMAGED, DO NOT PUT THE LENS BACK ON THE EYE. Place the lens in the storage case and contact the eye care practitioner.
- If the lens has dirt, an eyelash, or other foreign objects on it, or the problem stops and the lens appears undamaged, then thoroughly clean, rinse, and disinfect the lens; then reinsert it.
- After reinsertion, if the problem continues, the patient should IMMEDIATELY REMOVE THE LENS AND CONSULT THE EYE CARE PRACTITIONER.

When any of the above problems occur, a serious condition such as infection, corneal ulcer, neovascularization, or iritis may be present. The patient should be instructed to KEEP THE LENS OFF THE EYE AND SEEK IMMEDIATE PROFESSIONAL IDENTIFICATION OF THE PROBLEM AND PROMPT TREATMENT TO AVOID SERIOUS EYE DAMAGE.

**FITTING:**

Refer to the Professional Fitting Guide for detailed information on the fitting of the HDS Hybrid Contact Lens for daily wear. Copies are available from: SynergEyes, Inc. [www.synergeyes.com](http://www.synergeyes.com), Telephone: (USA) +1-760-476-9410

**RECOMMENDED INITIAL WEARING SCHEDULE:**

Although many practitioners have developed their own initial wearing schedules, the following sequence is recommended as a guideline. Patients should be cautioned to limit the wearing schedule to the level recommended by the eye care practitioner regardless of how comfortable the lenses feel. An initial daily wear schedule may be offered at the practitioner’s discretion; see example below:

Day 1: Wear not to exceed 6 hours total	Day 2: 6 hours
Day 3 - Day 5: 8 hours	Day 6: Wear as eye care practitioner allows during waking hours

The lens only may be worn for daily wear use and for a period of up to six (6) months or as recommended by the eye care practitioner. A well fit lens provides for centration and minimal movement. The effects of lid interaction on blinking and gravity may result in lens decentration during open eye wear. Studies have not been completed to show that the HDS Hybrid Contact Lens is safe to wear during sleep.

**Patients must be cautioned; “when in doubt, take it out.”** It is important that the new wearer not over wear the lens or endure a lens that has an obvious foreign body sensation. In the event of foreign body sensation, the patient should be instructed to remove the lens, clean and rinse it, and replace the lens on the eye. If the sensation continues, the lens should not be worn. The patient should report for follow-up evaluation at the prescribed follow up schedule. The visit is best scheduled after several hours of wear and the patient should report with the lens in place. This visit provides an excellent opportunity to evaluate lens fit, comfort, and vision. Upon the absence of clinical signs and complications, the patient may be instructed to continue daily wear of the lens, as instructed by the practitioner, until the next scheduled follow-up visit.

**WEARING AND LENS REPLACEMENT SCHEDULES SHOULD BE DETERMINED BY THE EYE CARE PRACTITIONER.**

**LENS CARE DIRECTIONS:**

Eye care practitioners should review lens care directions with the patient, including both basic lens care information and specific instructions on the lens care regimen recommended for the patient:

**General Lens Care (Clean and rinse, then disinfect lenses) – Basic Instructions:**

- Always wash, rinse, and dry hands before handling contact lenses.
- Always use fresh, unexpired lens care solutions.
- Use the recommended system of lens care, which is chemical (not heat), and carefully follow instructions on solution labeling. Different solutions cannot always be used together, and not all solutions are safe for use with all lenses. Do not alternate or mix lens care systems unless indicated on solution labeling.
- Do not use saliva, tap water, or anything other than the recommended solutions (e.g., non-preserved saline) for lubricating or rewetting lenses. Do not put lenses in the mouth.
- Lenses should be cleaned, rinsed, and disinfected each time they are removed. Cleaning and rinsing are necessary to remove mucus and film from the lens surface. Disinfecting is necessary to destroy harmful germs.

- Always remove, clean, rinse, enzyme (as recommended by the eye care practitioner), and disinfect lenses according to the schedule prescribed by the eye care practitioner. The use of an enzyme or any cleaning solution does not substitute for disinfection.
- Ensure that the tamper-proof seal on the solution container is intact prior to initial use; do not use if tamper-proof seal is broken or missing.
- Lenses should never be worn while swimming.
- Do not store the lens for prolonged periods at temperatures below 15° C or above 30° C.

NOTE: **Soft contact lens care products are recommended by SynergEyes for use with the SynergEyes HDS Hybrid Contact Lenses.** SynergEyes HDS (paflucofen D hem-iberfilcon A) Hybrid Contact Lenses **cannot** be cleaned and disinfected using the following contact lens solutions: OPTIFREE® PUREMOIST® Solution, Revitalens Ocutec® Multipurpose Disinfecting Solution, or Biotrue® Multipurpose Solution.

It is the responsibility of the practitioner to ensure that diagnostic lenses are properly cleaned, disinfected, rinsed, and stored between uses. Following the solution manufacturer’s recommendations for disinfection constitutes the minimum level of care necessary for disinfection as set by the FDA.

<b>1. Clean:</b>
While the recommended care products may be approved for a “No Rub” regimen, it is recommended that moderate daily cleaning be conducted for the HDS Hybrid Contact Lenses. Clean one lens first (always start with the same lens first to avoid mix-ups). Place the lens, front side down, in the palm of the hand and apply several drops of the contact lens cleaning or disinfecting solution. Using the ring finger of the other hand, apply slight pressure in a swirling motion for the time recommended by the disinfecting solution manufacturer. <b>Note: Do not clean the lens by rubbing it between the thumb and index fingers, as this may cause lens warpage.</b>

<b>2. Rinse:</b>
Rinse the lens thoroughly with the disinfecting solution to remove mucus and film from the lens surface; place the lens into the correct chamber of the lens storage case. Repeat this rinsing procedure for the second lens. Avoid rinsing the lens over the sink and avoiding the use of tap water as tap water contains many impurities that can contaminate or damage the lens and lead to eye infection or injury.

<b>3. Disinfect:</b>
After cleaning and rinsing the lenses, disinfect them by using the disinfecting solution. Follow the instructions provided with the disinfecting solution labeling. <b>NOTE: SynergEyes HDS Hybrid Contact Lenses cannot be heat (thermally) disinfected.</b> When using hydrogen peroxide lens care systems, instruct the patient to use ONLY the lens case provided with the hydrogen peroxide care system; this case is specially designed to neutralize the solution. Failure to use the specialized case will result in severe stinging, burning, and injury to the eye. Follow the recommendations on the hydrogen peroxide system labeling exclusively. Following disinfection with a peroxide system, the lenses should be rinsed with sterile saline prior to use. CAUTION: Lenses that are chemically disinfected may absorb ingredients from the disinfecting solution which may be irritating to the eyes. A thorough rinse in fresh, sterile rinsing solution prior to placement on the eye should reduce the potential for irritation.

<b>4. Store:</b>
To store lenses after disinfection, leave them in the closed/unopened case until ready to wear. Always keep lenses completely immersed in the disinfecting solution when the lenses are not being worn. If the patient discontinues wearing lenses but plans to begin wearing them again after a few weeks, they should ask the eye care practitioner for storage recommendations.

<b>5. Care of the Lens Case:</b>
Contact lens cases can be a source of bacteria growth. After removing the lenses from the lens case, empty and rinse the lens storage case with solution(s) recommended by the lens case manufacturer; then allow the lens case to air dry. When the case is used again, refill it with fresh storage solution. Lens cases should be replaced at regular intervals as recommended by the lens case manufacturer or the eye care practitioner.

<b>6. Lubricating/Rewetting Solutions:</b>
The eye care practitioner may recommend lubricating/rewetting solutions which can be used to wet (lubricate) lenses while they are being worn to make them more comfortable.

<b>7. Lens Deposits and Use of Enzymatic Cleaning Procedure:</b>
Enzyme cleaning may be recommended by the eye care practitioner. Enzyme cleaning removes protein deposits on the lens. These deposits cannot be removed with regular cleaners. Removing protein deposits is important for the well-being of the patient’s lenses and eyes. If these deposits are not removed, they can damage the lenses and cause irritation. <u>Enzyme cleaning does NOT replace routine cleaning and disinfecting.</u> Patients should carefully follow the instructions in the enzymatic cleaning labeling.

<b>8. Care for a Sticking Lens:</b>
If the lens sticks (stops moving) or cannot be removed, apply 5 drops of the recommended lubricating or rewetting solution directly to the eye and wait until the lens begins to move freely on the eye before removing it. To initiate movement gently push the lens upward with the index finger on the margin of the lower lid. If non-movement of the lens continues after 30 minutes, the patient should IMMEDIATELY consult the eye care practitioner.

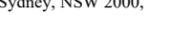
<b>9. Care for a Dehydrated Lens:</b>
The soft poly-hema portion of the HDS Hybrid Contact Lens may become dried out if left exposed to air while the lenses are off the eye. Patients should rehydrate the lens by carefully placing the lens into the storage case and covering it with the multipurpose solution. The lenses should be soaked for a minimum of five minutes prior to handling. Properly clean, rinse, and disinfect the lenses prior to reinsertion. If, after soaking, the lens does not become soft or if the surface remains dry, instruct the patient to NOT USE THE LENS UNLESS IT HAS BEEN EXAMINED BY THE EYE CARE PRACTITIONER.

**EMERGENCIES:**

If chemicals of any kind (e.g., household products, gardening solutions, laboratory chemicals, etc.) are splashed into the eyes, patients should: FLUSH EYES IMMEDIATELY WITH TAP WATER, THEN REMOVE LENSES PROMPTLY, IF POSSIBLE, AND IMMEDIATELY CONTACT THE EYE CARE PRACTITIONER OR VISIT A HOSPITAL EMERGENCY ROOM WITHOUT DELAY.

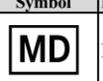
**REPORTING OF ADVERSE REACTIONS:**

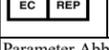
All serious adverse experiences and reactions observed in patients wearing HDS Hybrid Contact Lenses or experienced with the lenses should be reported to SynergEyes, Inc. and, as appropriate, the EU member state competent authority.

	SynergEyes, Inc., 2232 Rutherford Road Carlsbad, CA 92008, USA <a href="http://www.synergeyes.com">www.synergeyes.com</a> Phone: (760) 476-9410; Fax: (760) 476-9340	
	Emergo Europe: Westervoortsedijk 60, 6827 AT Arnhem, The Netherlands	
Australian Sponsor:	Emergo Asia-Pacific Pty. Ltd.: Level 20, Tower II, Darling Park, 201 Sussex Street, Sydney, NSW 2000, Australia	

**HOW SUPPLIED:**

Each HDS Hybrid Contact Lens is supplied sterile, in a sealed glass vial containing 0.9% buffered sodium chloride solution. NOTE: Do not dispense the lens if the tamper proof seal is broken or missing. The lenses are shipped as either a single lens or in a 2-pack carton package. The vial label is marked with the base curve radius/vault, skirt curve, diopter power, diameter, lot number, and expiration date. In addition, the vial label will show Rx Only and Sterile symbols. The packing slip or invoice is marked with the base curve radius, skirt curve, diopter power, diameter, and lot number.

Symbol	Definition	Symbol	Definition	Symbol	Definition
	Manufacturer		Country and Date of Manufacture		Medical device

Symbol	Definition	Symbol	Definition	Symbol	Definition
	Sterilized by irradiation; single sterile barrier system with protective packaging outside		Caution, consult instructions for use for important cautionary information		Caution: Federal law restricts this device to sale by or on the order of a licensed eye care practitioner.
	Lot number		Consult instructions for use		Unique device identifier
	Catalogue number		Use-by date		CE marking of conformity
	Upper/Lower limit of temperature		Avoid contact with water		Authorized representative in the European Community
	Do not use if package is damaged		Single patient multiple use	<u>Parameter Abbreviations:</u> Diameter (DIA), Base Curve Radius (BC), Vault (VLT), Skirt Curve Radius (SC), Diopter Power (PWR), Add Power – Multifocal Only (ADD), Zone – Multifocal Only (ZN), Lift – PS Lens Only (LIFT), Center Thickness (CT)	

**SUMMARY OF CLINICAL STUDIES:**

**SynergEyes A and SynergEyes M:** A three-month clinical study of the SynergEyes Hybrid Contact Lens was conducted to assess safety and effectiveness for vision correction in daily wear that included subjects with nearsightedness, farsightedness, and astigmatism. Included in this population were 13 subjects (13 eyes) that fitted for monovision in one eye for near vision correction. The study was designed to evaluate contact lens visual acuity and wearing time, and assess contact lens adverse events and loss of visual acuity. Overall Findings: A total of 107 subjects in which 61 (57%) subjects completed, and 46 (43%) subjects were discontinued. The population demographics were similar to previous contact lens studies with a female to male gender ratio of 2.3 to 1.0. The average age of the completed and discontinued subjects was 42.1 with an age range of 12 years to 60 years of age. Safety: There were no serious adverse events reported during the trial. There was one subject (2 eyes) discontinued for complications of corneal abrasions, described as “poor outcome with lenses”. The lens was treated, and the condition resolved in both eyes. Additional reasons for discontinuation included poor comfort (32.6%), loss of interest (28.3%), poor vision (13.0%), handling difficulty (2.2%), and other (23.8%). “Other” reasons included unrelated medical conditions (4.3%), non-compliance (8.7%), poor outcome with lenses (4.3%) including the above-described subject, and undefined (6.5%). Efficacy: Visual Acuity- Final visual acuity for completed subjects was 20/20 or better (40.2%), 20/25 or better (65.6%), 20/30 or better (82.0%), and 20/40 or better (90.2%). No visual acuity was lost during the course of the study. Visual acuity was reported worse than 20/40 for 9.8% (12/122) of the completed eyes. Of those 8 eyes were fitted for near vision for monovision correction. Wearing Time: Over the study period, the average daily wearing time reported by completed patients was 11.4 hours per day. Conclusion: The SynergEyes Hybrid Contact Lens provided satisfactory performance as expected. Discontinuation rates were somewhat higher than normal due to the reasons cited above. There were no discontinuations for safety related reasons. Overall, the lens performance demonstrated safe and effective use of the device for its intended use.

**SynergEyes KC and ClearKone:** A one-month clinical study of the SynergEyes KC Hybrid Contact Lens was conducted to assess safety and effectiveness for vision correction in daily wear for patients suffering from keratoconus and associated refractive errors of nearsightedness, farsightedness, and irregular astigmatism. The study was designed to evaluate contact lens visual acuity and wearing time; and assess contact lens adverse events and loss of visual acuity. Overall Findings: A total of 62 subjects were dispensed into the study of which 44 (71%) subjects completed, and 18 (29%) subjects were discontinued. The population demographics were similar to other contact lens studies with a female to male gender ratio of 1.48 to 1.0. The average age of the completed participants was 40.7 and the average age of the discontinued participants was 39.4. Safety: Four (4) adverse events were reported during the study for 2 completed subjects and 2 discontinued subjects. One subject presented with mechanical abrasion on the apex of the cone in both eyes after several lens dispensing visits, one subject experienced transient changes in intraocular pressures in both eyes at the 1-week follow-up visit, one subject presented with edema and infiltrates in the left eye 61 days after first lens dispensing, and one subject presented with a red left eye with small infiltrates 14 days after first lens dispensing. Eighteen (18) subjects, or 29% of the 62 dispensed subjects, were discontinued. The most common reasons for discontinuation were poor comfort (38.9%), poor outcome with lenses (16.7%), and ‘other’ (16.7%). “Other” reasons included patients moved, and lens ripping- too time consuming to restart. Additional reasons for discontinuation included poor vision (11.1%), non-compliance (11.1%), and loss of interest (5.6%). The most common symptoms, problems, and complaints for completed eyes were cited as discomfort/awareness (26%), blurred vision (11.7%) and dryness-scratchiness (12.5%). These rates were higher for discontinued patients (41.6%, 24.7%, and 18.7% respectively). Efficacy: Visual Acuity- Final visual acuity for completed subjects was 20/20 or better (26.3%), 20/25 or better (55.1%), 20/30 (71.4%), and 20/40 or better (83.9%). The visual acuity rates for discontinued subjects were 8.4%, 27.8%, 47.2%, and 52.8%, respectively. Vision correction fluctuated as expected with the instability of the corneal curvature from keratoconus under the contact lens contributing to the change. Five (5) completed eyes and 3 discontinued eyes were reported to have VA decreases of more than 2 lines of Snellen VA when comparing the contact lens VA with the best corrected VA. These findings are expected with this population. Conclusion: The *SynergEyes KC* Hybrid Contact Lens for Keratoconus provided satisfactory performance as expected. Discontinuation rates were somewhat higher than normal due to the nature of the subjects enrolled in the study. There were no discontinuations for safety related reasons. Overall, the lens performance demonstrated safe and effective use of the device for Keratoconus.

**SynergEyes PS:** A three-month clinical study of the SynergEyes PS Hybrid Contact Lens for Post-Surgical Refractive Error and Trauma was conducted to assess safety and effectiveness for vision correction in daily wear for patients suffering from post-surgical refractive error including nearsightedness, farsightedness, astigmatism and irregular astigmatism, or trauma to the eye. The study was designed to evaluate contact lens visual acuity and wearing time; and assess contact lens adverse events and loss of visual acuity. Overall Findings: A total of 76 subjects were dispensed into the study of which 44 (80 eyes) [57.9%] completed the study and 32 subjects (58 eyes) [42.1%] were discontinued from the study. An additional 8 subjects were not dispensed lenses. The population demographics were similar to other contact lens studies with a female to male gender ratio of 1.3 to 1.0. The average age of the completed and discontinued subjects was 49 and 43.5, respectively. Safety: Eight (8) adverse events were reported during the study for 4 completed subjects and 3 discontinued subjects. The adverse events included: 1 painful light sensitivity (photophobia); 1 loose interrupted suture and infiltrate; 1 abrasion of the cornea upon lens removal; 1 small superficial abrasion; 1 subject (2 eyes) with swelling of the corneas due to inadequate rinsing of hydrogen peroxide disinfection solution from lenses; 1 allergic conjunctivitis; and 1 keratitis iritis. Thirty-two (32) subjects (58 eyes) [42.1%] were discontinued from the study. An additional 8 subjects were never dispensed lenses. The most common reasons for discontinuation were poor outcome with lenses (32.5%), poor comfort (27.5%), poor vision (12.5%), loss of interest (10.0%), non-compliance (10.0%), and handling difficulty (2.5%). Symptoms, Problems, and Complaints: For the completed eyes, no symptoms were reported at 46.2% of the dispensing visit or follow-up visit examinations and for discontinued eyes no symptoms were reported at 28.8% of the dispensing visit or follow-up visit examinations. Symptoms decreased over time for completed eyes but not for discontinued eyes. The most common symptoms reported were discomfort and awareness (28.4% completed, 45.8% discontinued), dryness and scratchiness (23.4%, 22.1%), itchiness and burning (9.3%, 8.9%), and variable vision (8.2%, 13.3%). Efficacy: Visual Acuity- Final visual acuity for completed subjects was 20/20 or better (20.0%), 20/25 or better (45.0%), 20/30 (68.8%), and 20/40 or better (81.3%). The visual acuity rates for discontinued subjects were 8.6%, 39.7%, 53.4%, and 70.7%, respectively. Vision correction fluctuated as expected with the instability of the corneal curvature from post-surgical conditions and trauma under the contact lens contributing to the change. Three (3) completed eyes and 8 discontinued eyes were reported to have VA decreases of more than 2 lines of Snellen VA when comparing the contact lens VA with the best corrected VA. These findings are expected with this population. Wearing Time: Over the study period the average daily wearing time reported by completed patients was 10.6 hours per day. Conclusion: SynergEyes PS Hybrid Contact Lenses for post-surgical refractive error and cases involving trauma provided satisfactory performance as expected. The higher than estimated discontinuation rate was anticipated due to the nature of the subject population with the inclusion of subjects who might otherwise have less successful outcomes with other lens types. Overall, the lens performance demonstrated safe and effective use of the device for its intended use.